


INSTRUCTIONS FOR USE

AbbVie Cross Body Pack


For use with CADD-Legacy® 1400 pump

This guide pertains only to the use of the AbbVie Cross Body Pack. For questions or problems, refer to these Instructions for Use or call support toll free at 1-844-386-4968.

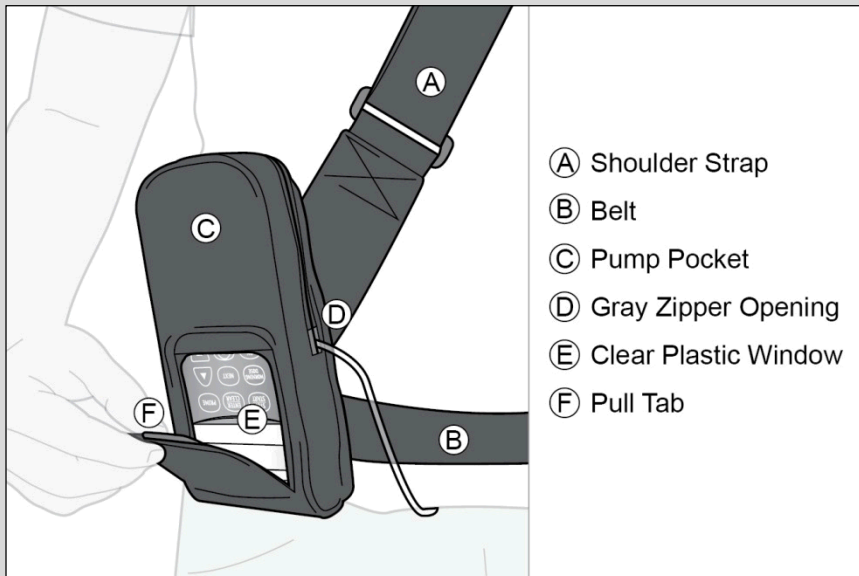

About your Cross Body Pack:

The Cross Body Pack is designed to carry a CADD-Legacy® 1400 pump. This Cross Body Pack is intended for single-patient use only. The Cross Body Pack should be worn over the shoulder and positioned at or near the waist using the Belt for support (see Figure i).

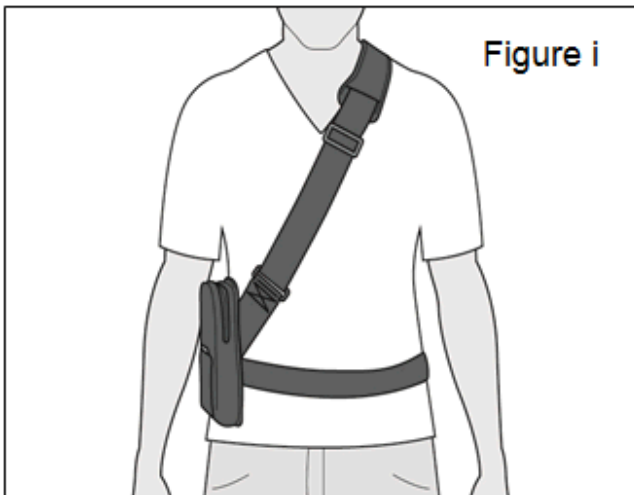


Caring for your Cross Body Pack:

Care: Do not machine wash.

Precautions:

- ⚠ **CAUTION:** When putting on the Cross Body Pack, ensure the pump is resting on a flat surface, and will not fall (see Figure ii). **If the pump is not completely secure, it may fall and injure you or damage the pump.**

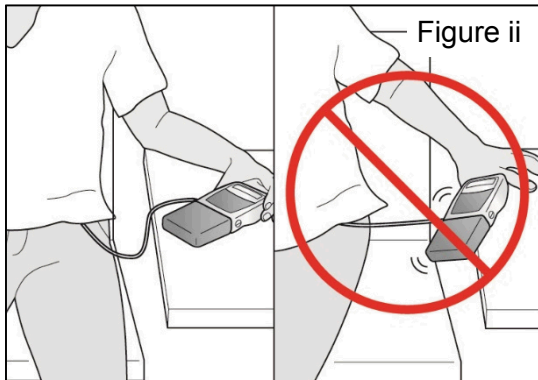


- ⚠ **CAUTION:** In the event the pump is dropped and no noticeable injury is suffered, inspect the pump for damage. In the event of pump damage contact your healthcare provider.

To begin using your Cross Body Pack:


1) Put on the Cross Body Pack:

⚠ **CAUTION:** Make sure that the pump is stable on a flat surface, and will not fall while you put on the Cross Body Pack.

a. Hold the Cross Body Pack so that the Pull Tab faces out and the zipper is on top (see Figure A).

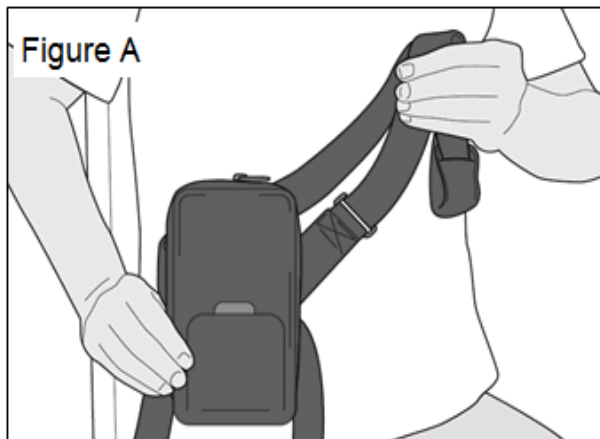


b. Put your arm through the Shoulder Strap and then place the Shoulder Strap over your head (see Figure B).

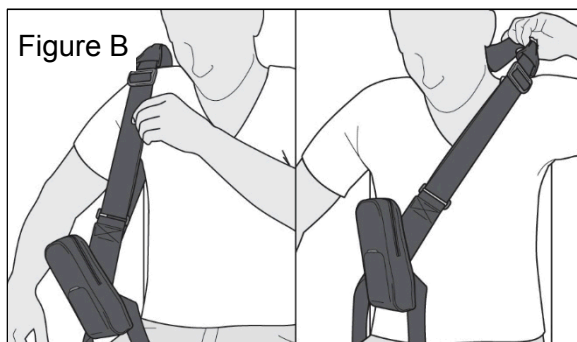

c. Confirm that the Shoulder Strap lies flat against your body and the Shoulder Pad rests comfortably on your shoulder.


2) Put on the Belt:

- a. Wrap the Belt around your waist and confirm that the Belt lies flat against your body.
- b. Fasten the clip. Listen for a snap sound.
- c. Check that the clip is securely fastened (see Figure C).

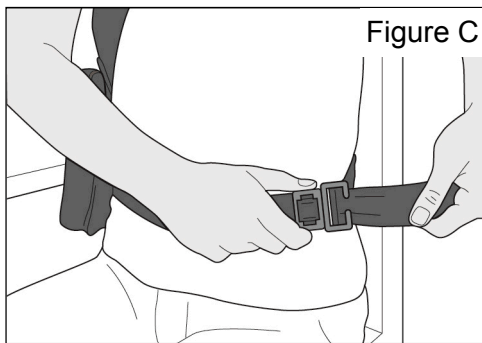

3) Adjust the Cross Body Pack:


- a. Move the plastic slider until the Shoulder Strap is at the desired length (see Figure D).

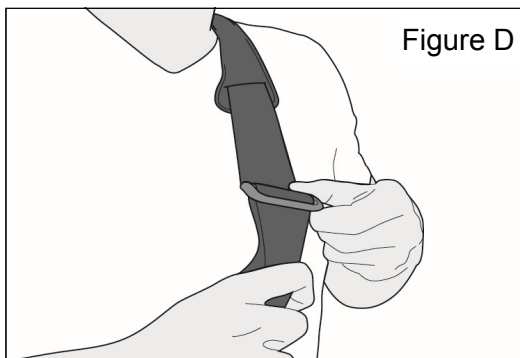

- b. Move the plastic slider until the Belt is at the desired length. The Cross Body Pack should fit comfortably but not move or shift.

4) Insert the pump into the Pump Pocket:


- a. Place the pump into the Pump Pocket so that the screen is on the bottom and facing out (see Figure E).

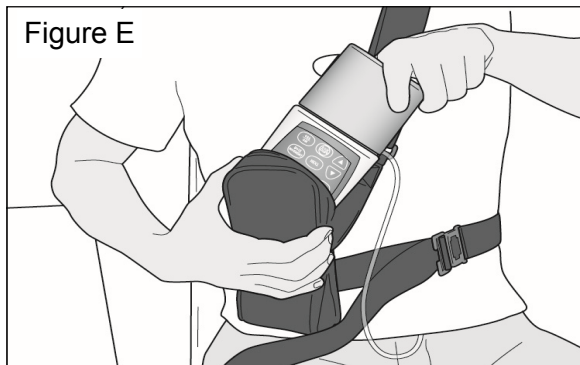

- b. Using the Pull Tab, confirm that both the display and the keypad are visible through the Clear Plastic Window.
- c. Route the tubing to exit through the Gray Zipper Opening on the side of the Pump Pocket. Make sure that the tubing is not kinked.
- d. Zip the Pump Pocket shut.

5) Interact with the pump (when needed):

- a. You can access the pump's keypad through the Clear Plastic Window.
- b. Use the Pull Tab to view or cover the Clear Plastic Window (see Figure F).

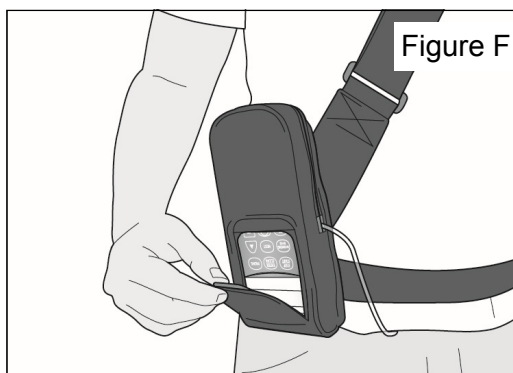

To take off your Cross Body Pack:

1) Remove the pump from the Pump Pocket:

⚠ CAUTION: Make sure that the pump is stable on a flat surface, and will not fall while you take off the Cross Body Pack.

- a. Pull the zipper to open the Pump Pocket.
- b. Remove the pump from the Pump Pocket.
- c. Place the pump on a flat surface.

2) Take off the Cross Body Pack:

- a. Undo the clip on the Belt.
- b. Remove the Shoulder Strap over your head.

Disposing of your Cross Body Pack:

Dispose of the device according to your local government ordinances and recycling plans. Otherwise contact AbbVie for help with device disposal. Each device is for single patient use only and cannot be re-sold or used on multiple patients.

Explanation of Symbols (per ISO 15223-1):



Catalog number



Batch code



Manufacturer



Medical Device



Use-by date



rxabbvie.com

Consult *instructions for use* or consult *electronic instructions for use*



Caution



10636-04-01



Manufactured for:
AbbVie Inc.
1 North Waukegan Road
North Chicago, IL 60064 USA



rxabbvie.com

rxabbvie.com

©2021 AbbVie Inc. All rights reserved.